
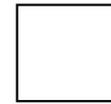


ENGINEERING SPECIFICATIONS: Vent shall be manufactured from polypropylene resin and shall conform to ASTM #1412-92 for corrosion resistance and material content. Vent will contain no PVC or solvents used for assembly. Vent may be constructed using rigid or flexible polypropylene pipe as required. Vent will be manufactured by Scientific Plastics Company of Kansas City, KS, USA.

NOTES:


- Assembly reference: Compression fittings are to be handtight only.
- Flex tubing is to be installed with an upward draft with no kinks or loops for maximum airflow.
- Rigid pipe and flex tubing is easily trimmed to fit in the field.
- Minimum space required to accommodate sweep 4.55 from cabinet wall.
- All dimensions are in inches.
- Custom vents & lengths are available.

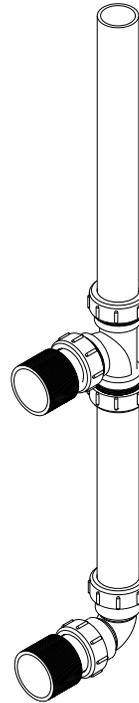

Check box for optional squeezed Top penetration. For oval cutout With dims of 1.965 X 1.437.



Oval Cutout

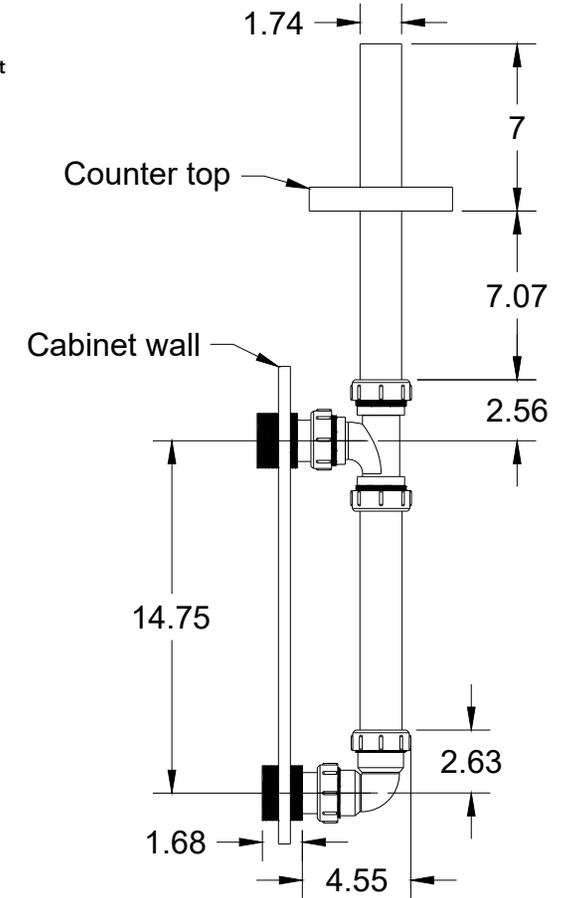

Narrow chase vent



(Installed with S.P. cabinet liner)

MODEL: F160L089-150		PRODUCT SPECIFICATION	
DESCRIPTION: Double solvent rigid vent - standard top penetration - narrow chase.		SCIENTIFIC PLASTICS COMPANY	
DWG BY: Philip Llamas		DATE: 12/19/2018	
		 1016 SOUTHWEST BLVD KANSAS CITY KS 66103 (913)-432-0322 OFFICE (913)-432-0459 FAX 	